


# FORM FIT®

## Hydrogel Intracanalicular Plug

Form Fit® comes sterile, preloaded in an inserter. The polyimide sheath is inserted into the punctum and the inserter handle compressed to expel the hydrogel plug. Once inserted the plug increases in size as it hydrates until it completely fills the vertical canalicular cavity to form a perfect fit.

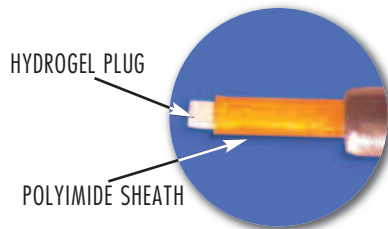


### PRODUCT FEATURES


- \* In its dry state, the plug is 0.3mm in diameter by 3.0 mm long and requires no special storage
- \* In the hydrated state, the plug is a soft, pliable, gelatinous material
- \* The need for punctal measurement has been eliminated
- \* Hydrogel is biocompatible and causes no foreign body sensation

REF #	Description	Qty per Box
6303	Form Fit® Hydrogel Plugs	10 plugs
6303-T	Form Fit® Hydrogel Plugs	2 plugs

INSERTION

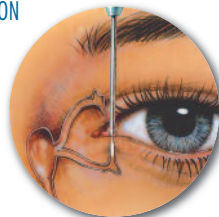

IMMEDIATELY AFTER INSERTION

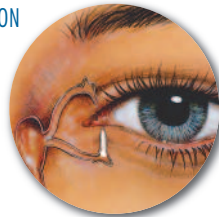

10 MINUTES AFTER INSERTION

